

(For scientific research use only, not for clinical diagnosis!)

Monkey Hyaluronic Acid (HA)

ELISA Kit Instructions for Use

Specifications: 48T/96T

Please read the instructions carefully before use. If you have any questions,

please contact us through the following methods: Official hotline: 025-5229-

8998 Sales department phone: 13914481711 Technical phone: 15950492658

Contact email: 3224949330@qq.com Company website:

www.byabscience.cn For specific shelf life, please see the reagents Box

packaging label. Please use the kit within the shelf life.

When contacting us, please provide the produ	ct number and production	date (see box	label) so that w	_{re}
can serve you more efficiently.				

Nanjing BYabscience technology Co.,Ltd



Kit performance Detection range: 0.625

ng/ml-10 ng/ml.

Sensitivity: The lowest detectable dose is less than 0.1 ng/ml.

Precision: intra-batch variation coefficient CV% is less than 10%; inter-batch variation coefficient CV% is less than 15%.

Recovery rate: The recovery rate is between 85%-115%.

Specificity: This kit recognizes natural and recombinant monkey hyaluronic acid (HA) and has no crossover with structural analogs.

Stability: Stored at 2℃-8℃, validity period is 6 months.

Purpose: Used to detect the concentration of monkey hyaluronic acid (HA) in samples such as serum, plasma, cell culture supernatant and tissue. Shelf life: Stored at 2°C-8°C, valid for 6 months.

Experimental principle

The kit uses enzyme-linked immunoassay method. Biotin-labeled HA is used, and the purified anti-HA antibody is coated on the microplate. In the competitive inhibition reaction, a certain amount of solid-phase antibody is used to inhibit the competitive reaction with biotin-labeled HA and non-labeled antigen (calibrator or specimen). The antibody The amount of binding to biotin-labeled HA is inhibited by the amount of non-labeled antigen. The greater the amount of non-labeled antigen, the less the antibody will bind to biotin-labeled HA, and vice versa. After the reaction is balanced, a solid-phase antibody-biotinylation is formed. HA, and then add enzyme-labeled avidin to form a solid-phase antibody-biotinylated HA-enzyme-labeled-avidin complex. After adding substrate for color development, use a microplate reader to measure the absorbance (OD value) at a wavelength of 450 nm. As the HA concentration increases, the OD value gradually decreases with a good linear

relationship. This kit has the characteristics of high sensitivity, strong specificity, good repeatability, simple and rapid operation, and has reliable detection performance for the reduction or increase of HA in serum.

Nanjing BYabscience technology Co.,Ltd

Kit components and storage: Store unopened kits at 2-8

degrees Celsius. Do not use expired kits.

Components	48-well configuration	96-well configuration	Store after opening
Pre-coated enzyme	48T	96T	2-8°C14 days
Standard product	0.3mL*6 tubes	0.3mL*6 tubes	2-8°C14 days
sample diluent	3ml	6ml	2-8°C180 days
biotinylated antigen	3ml	6ml	2-8°C14 days
HRP labeled avidin	3ml	6ml	2-8°C180 days
Chromogenic substrate	3ml	6ml	2-8°C180 days
Chromogenic substrate	3ml	6ml	2-8°C180 days
stop solution	3ml	6ml	2-8°C180 days
20×Lotion	15ml	25ml	2-8°C180 days
sealing film	2 sheets	2 sheets	
manual	1 serving	1 serving	
Ziplock bag	1	1	

The concentrations of calibrators are: 10, 5, 2.5, 1.25, 0.625, 0 ng/ml.

Note: 1: Before use, please check whether the label and quantity of the reagents in the kit are consistent with the table.

2: If the components of the kit need to be used again, please ensure that they have not been contaminated since the last use. 3: If the enzyme plate is not used up in a single time, remember to seal it and store it at 2-8°C.

Prepare your own test equipment required for the test (not provided, but can assist in

1) Microplate reader capable of detecting absorbance at 450 nm 2) Pipette, pipette tip, and sample addition tank 3) 37°C incubator or water bath 4) Test tubes, centrifuge tubes, measuring cylinders, etc. for preparing reagents 5)

Distilled water or deionized water Ionized water

Nanjing BYabscience technology Co.,Ltd



6) Vortex oscillator and microplate oscillator.

Notes: 1) For scientific research use only, not for clinical diagnosis.

- 2) Use within the validity period marked on the kit. Expired products must not be used.
- 3) Do not mix with kits or components from other manufacturers. Use the sample diluent provided with the kit.
- 4) If the sample value is higher than the highest standard concentration value, please dilute the sample appropriately and then re-measure.
- 5) Human anti-mouse and other heterophilic antibodies present in the sample to be tested will interfere with the test results. Please eliminate this factor before testing.
- 6) The test results obtained by other methods are not directly comparable to the test results of this kit.
- 7) Please wear a lab coat and latex gloves for protection during the test. Especially when testing blood or other body fluid samples, please follow the national biological laboratory safety protection regulations.
- 8) Carry out incubation strictly according to the specified time and temperature to ensure accurate results. All reagents must reach room temperature 20-25°C before use. Store reagents refrigerated immediately after use.
- 9) Improper plate washing can lead to inaccurate results. Make sure to absorb as much liquid as possible from the wells before adding substrate. Do not allow the microwells to dry out during incubation.
- 10) Eliminate residual liquid and fingerprints on the bottom of the plate, otherwise it will affect the OD value.
- 11) The substrate chromogenic solution should be colorless or very light in color.
- 12) Avoid cross-contamination of reagents and specimens to avoid erroneous results.

- 13) Avoid direct exposure to strong light during storage and incubation.
- 14) The microplate reader used for detection needs to be equipped with a filter capable of detecting a wavelength of 450±10nm, and the optical density range is between 0-3.5. It is recommended to preheat 15 minutes in advance before use.
- 15) The EP tubes and tips used in the test are single-use and are strictly prohibited from mixing.

Nanjing BYabscience technology Co.,Ltd



Sample preparation and storage

The following lists only general guidelines for sample collection and preservation. During the collection and storage of all samples, sodium azide must not be used as a preservative. If the sample is not analyzed immediately, it should be aliquoted and stored frozen, and repeated freezing and thawing should be avoided.

Cell culture supernatant - centrifuge to remove precipitate, analyze immediately or aliquot and store frozen at -20°C.

Serum - Collect blood in a clean test tube, coagulate at room temperature for 30 minutes, centrifuge at 2000×g for 20 minutes, and collect serum. Analyze immediately or aliquot and store frozen at -20°C.

Plasma—anticoagulate with heparin, citrate, or EDTA, and centrifuge at 2000×g for 20 minutes at 2-8°C within 30 minutes of blood draw. To eliminate the influence of platelets, it is recommended to further centrifuge at 10,000 × g for 10 minutes at 2-8°C. Analyze immediately or aliquot and store frozen at -20°C.

Cell lysis buffer - For adherent cells, remove the culture medium and wash with PBS, normal saline or serum-free culture medium. Add an appropriate amount of lysis solution and pipet several times with a gun to fully contact the lysate and cells. Typically after 10 seconds, cells are lysed. For suspended cells, collect the cells by centrifugation and wash them with PBS, physiological saline or serum-free culture medium. Add an appropriate amount of lysis solution, blow the cells with a gun, and flick them with your fingers to fully lyse the cells. After full lysis, centrifuge at 10000-14000×g for 3-5 minutes and take the supernatant. Analyze immediately or aliquot and store frozen at -20°C.

组织匀浆——用预冷的 PBS (0.01M, pH=7.4)冲洗组织,去除残留血液(匀浆中裂解的红细胞会影响测量结果),称重后将组织剪碎。将剪碎的组织与对应体积的 PBS (一般按 1:9 的重量体积比,比如 1g 的组织样品对应 9mL 的 PBS,具体体积可根据实验需要适当调整,并做好记

录。推荐在 PBS 中加入蛋白酶抑制剂)加入玻璃匀浆器中,于冰上充分研磨。为了进一步裂解 组织细胞,可以对匀浆液进行超声破碎,或反复冻融。最后将匀浆液于 5000×g 离心 5~10分钟,取上清检测。

尿液——用无菌管收集, 离心 2000×g 20 分钟。仔细收集上清。如有沉淀形成, 应再次离心。

Nanjing BY abscience technology Co., Ltd

网址: www.byabscience.cn 官方热线: 025-5229-8998 监督电话: 15950492658



试剂准备 1、使用前,所有的组分都要至少复温 60min,确保充分复温到室温。

2、浓缩洗涤液: 从冰箱取出的浓缩洗涤液, 会有结晶产生, 这属于正常现象, 水浴加热使结晶完全溶解。浓缩洗涤液与蒸馏水, 按 1:20 稀释, 即 1 份的浓缩洗涤液, 添加 19 份的蒸馏水。

操作程序所有试剂和组分都先恢复到室温,标准品、质控品和样品,建议做复孔。

- 1、按前面说明书描述的方法,配制好试剂盒各种组分的工作液。
- 2、从铝箔袋中取出所需板条,剩余的板条用自封袋密封放回冰箱。
- 3、将预包被板从密封袋中取出,设一个空白对照孔,不加任何液体;每个校准品设 2 孔,每孔加入对应校准品 50 ul;其余每个检测孔直接加待测血清或质控品 50 ul。
- 4、除空白孔外所有孔加入生物素化抗原 50μl,混匀,贴上封板膜,置 37℃温育 60 分钟。
- 5、手工洗板:弃去孔内液体,洗涤液注满各孔,静置 10 秒甩干,重复 3 次后拍干。洗板机洗板:选择洗涤 3 次程序洗板后拍干。

(提示: 为获得理想的实验结果,必须彻底移除残留液体。洗板完成之后,请立即进行下步操作,不要让微孔板干燥。)6、每孔加入酶标亲和素 50μl(空白对照孔除外),混匀,贴上封板膜,置 37℃温育 30 分钟。

- 7、手工洗板:弃去孔内液体,洗涤液注满各孔,静置 10 秒甩干,重复 3 次后拍干。洗板机洗板:选择洗涤 3 次程序洗板后拍干。
- 8、每孔加显色剂 A 50μl, 显色剂 B 50μl, 振荡混匀后,置 37℃避光显色 15 分钟,每孔加终止液 50μl。

9、用酶标仪读数,取波长 450nm, 先用空白对照孔调零点, 然后测定各孔光密度值(OD值)。


Nanjing BYabscience technology Co.,Ltd

网址: www.byabscience.cn 官方热线: 025-5229-8998 监督电话: 15950492658



结果计算

9、以标准品浓度做为横坐标,对应的吸光度(OD 值)作为纵坐标,利用计算机软件,采用四参数 Logistic 曲线拟合(4-pl),创建标准曲线方程,通过样本的吸光度(OD 值),利用方程计算样品的浓度值。【用 ELISA Calc 软件计算】10、如果样品被稀释,通过上述方法测的的浓度值,要乘以稀释倍数,才是样品的最终浓度。注意:实验者需根据自己的实验建立标准曲线。每次检测,每块酶标板都必须设立标准曲线。以下曲线仅供参考!

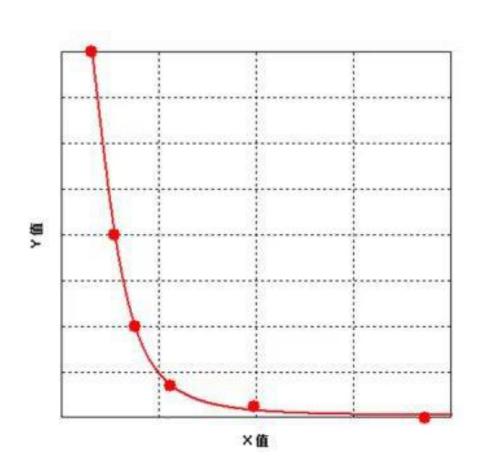

(标曲示意图,仅供参考)

Nanjing BYabscience technology Co.,Ltd

网址: www.byabscience.cn 官方热线: 025-5229-8998 监督电话: 15950492658



[问题分析] 若实验效果不好,请及时对显色结果拍照,保存实验数据,保留所用板条及 未使用试剂,然后联系我公司技术支持为您解决问题。同时您也可以参考以下资料: [问题解答]

问题描述	可能原因	相应对策相应对策	
	吸液或加液不准	检查移液器及吸头	
标准曲线梯度差	平衡时间太短	保证充足的平衡时间	
	洗涤不完全	保证洗涤时间和洗涤次数及每孔的 加液量	
显色很弱或无色	孵育时间太短	保证充足的孵育时间	
	实验温度不正确	使用推荐的实验温度	
	试剂体积不够或漏加	检查吸液及加液过程,保证所	
	稀释不正确	有试剂按顺 序足量添加	
	酶标记物失活或底物失效	混合酶结合物和底物,通过迅速显色来检查判断	
读数数值低	酶标仪设置不正确	在酶标仪上检查波长及滤光片设 提前打开酶标仪预热	
变异系数大	加液不正确	检查加液情况	
背景值高	检测抗体的工作浓度过高	使用推荐的稀释倍数	
	酶标板洗涤不完全	保证每步清洗完全;如果用自动洗板机,请检查所有的出口 是否有堵塞;是否使用试剂盒 配备的洗涤液	
	洗液有污染	配制新鲜的洗液	
灵敏度低	Improper storage of ELISA kits	Store relevant reagents according to	
7	Not terminated before reading	Stop solution should be added to	